

UNICRYL quick

5 days

*Effective Tensile
Strength*


UNIMED[®]
SUTURES

United Medical Industries Co. Ltd.

www.unimed.com.sa

UNICRYL quick

- Coated Braided Polyglycolic Acid B.P. Undyed and Violet.
- Synthetic rapidly absorbed sterile surgical suture composed of a homopolymer of glycolide.
- Available in range of gauge sizes and lengths, non-needed or attached to stainless steel needles of varying types and sizes.
- Minimal tissue reaction.
- Easy passage through tissue.
- Secure Knot tying and easy knot tie down.
- Intended for use where an absorbable suture or ligature is indicated for short term support. It can be used for superficial soft tissue approximation of the skin and mucosa.
- Sterilised by Ethylene Oxide.
- Coating – Polyglycolic acid is a homopolymer of glycolic acid, which is braided into sutures of various sizes. The braid is coated with mixture of Calcium Stearate and Polycaprolactone.
- Sutures are guaranteed to meet the requirements of the British Pharmacopoeia (B.P.), the European Pharmacopoeia (Ph.Eur.), and the United State Pharmacopoeia (U.S.P.).

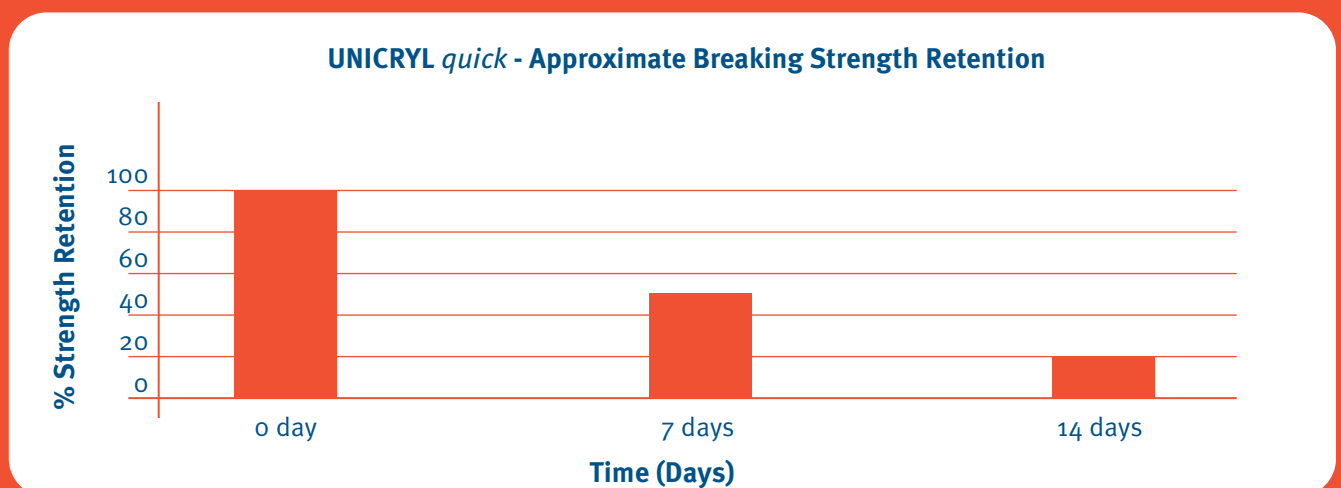

Head Office

P.O.Box 51871
Riyadh 11553 Saudi Arabia
Tel. +966 1 2650665
Fax. +966 1 2650569
general@unimed.com.sa